

Oxygen sensors in AquaResp

Pyroscience examples
by
Morten Bo

Pyroscience Firesting Firmware >4

For the latest firmware for Firestings, you need to use AquaOxyLog to communicate oxygen measurements to Aquaresp.


This is easily done by added the path of the AquaResp oxygen folder to OxyLog settings.

See next slide for images and arrows.


When settings are done. Start OxyLog, then Aquaresp.

OxyLog and AquaResp communication

AquaOxyLog Directory

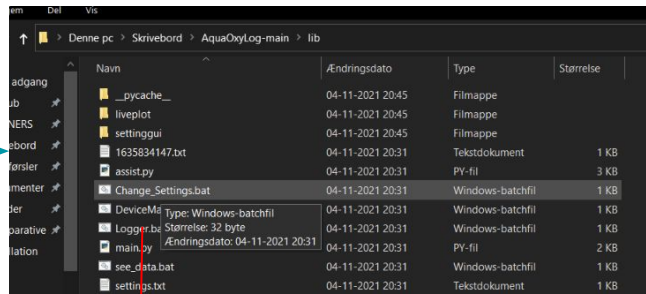

1 Open Change_Settings.bat
Can also be opened from the
settinggui folder > open HTML file

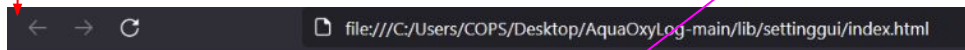

Create settings file for AquaOxyLog

This form **does not** show the present settings. It can only create new files


Active channels:

Location datafile:

Location aquaresp:



Device communication port:

Press save

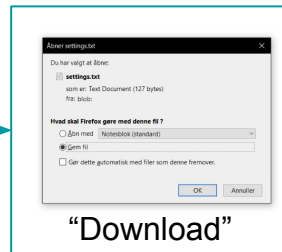

"Download"

AquaResp OXYGEN Directory

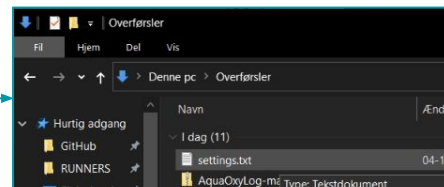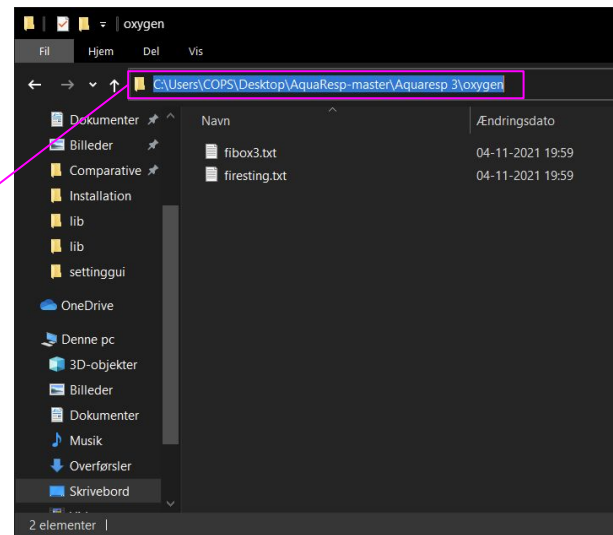

Move or copy Settings.txt to OxyLog Directory and REPLACE previous

Note


Active channels in OxyLog should match that of AquaResp

Settings can be change via the bat file, but also using the html file in ../settinggui
the settings file should be named “settings.txt” for OxyLog to find it

On starting OxyLog

Start Aqua OxyLog


1

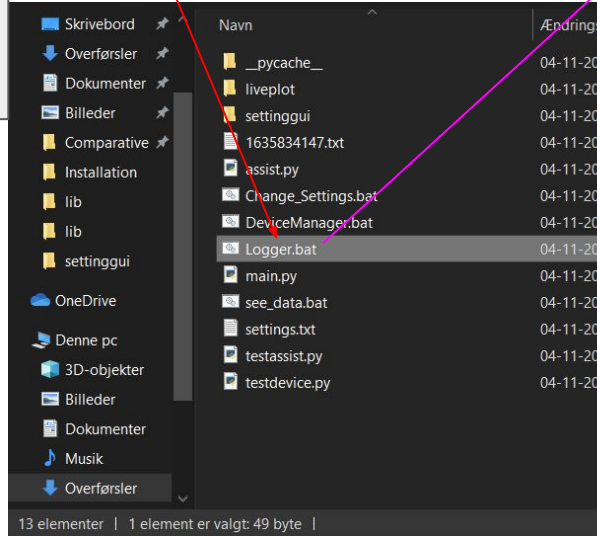

```
C:\WINDOWS\system32\cmd.exe
Python was not found; run without arguments to install from the Microsoft Store, or disable this shortcut from Settings
> Manage App Execution Allases.
C:\Users\COPS\Desktop\
Channel: 1 20.789 1636054888.244424
Channel: 2 21.178 1636054888.9102476
C:\Users\COPS\Desktop\
Channel: 1 20.796 1636054889.6665134
Channel: 2 21.188 1636054890.2653296
C:\Users\COPS\Desktop\
Channel: 1 20.803 1636054891.011121
Channel: 2 21.192 1636054891.6149619
C:\Users\COPS\Desktop\
```

2


Check AquaResp Oxygen folder

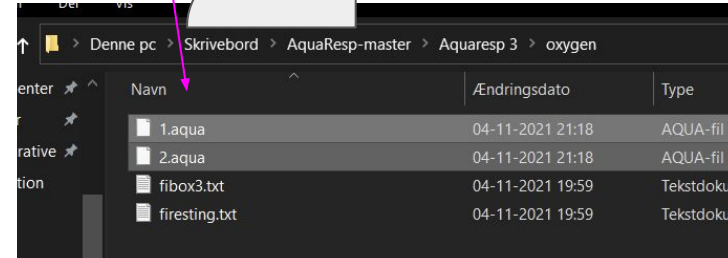



2 active channels should give 1.aqua and 2.aqua
The format is just a text file. Open in notepad to see contents if you are interested

Open AquaResp

Open AquaResp. The default choice for a new installation is Firesting New but go to experimental settings. Make sure **Firesting New** is chosen.

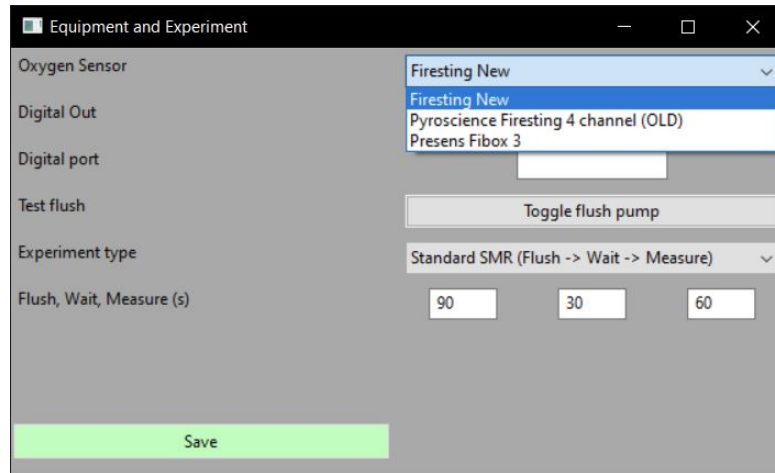

Check that measurements are live in AquaResp

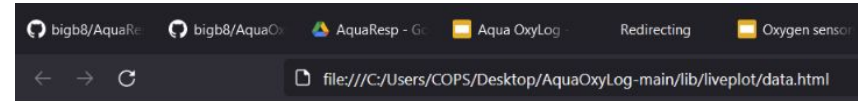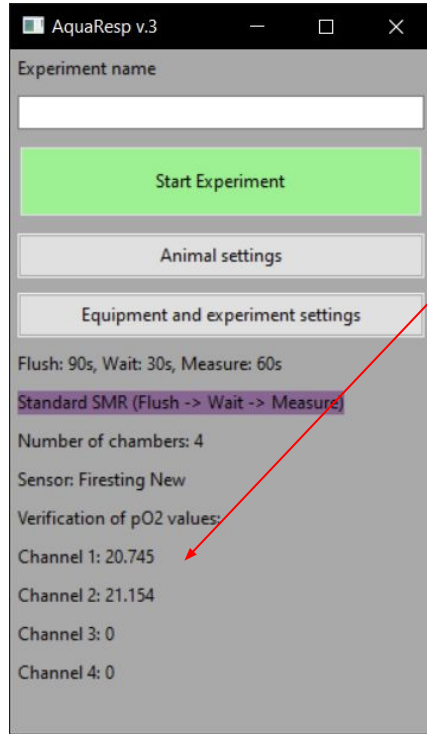

Channel 1

Oxygen: 20.73 [%air sat]
at 2021-11-04 20:24:58
(yyyy-mm-dd hh:mm:ss)

Channel 2

Oxygen: 21.143 [%air sat]
at 2021-11-04 20:24:58
(yyyy-mm-dd hh:mm:ss)

Measurements in the html gui on OxyLog
are refreshed every 5s. But sampled approx
1s

Firmware < 4 (per nov 2021)

1 Choose Firesting ... (OLD)

2 Open Firesting SW and make
it save to AquaResp
3/oxygen/firesting.txt

3 Check that AquaResp
measurements are live

